

## COMPOSITIONS AND METHODS FOR DETECTING RAPHIDOPHYTES

### RELATED APPLICATION

5 This application is a Continuation in Part Application of application serial  
number 09/596,136 filed June 16, 2000 which claims' priority under 35 U.S.C. §  
119(e) of United States provisional patent application number 60/141,362 filed  
June 28, 1999, each of which are hereby incorporated by reference in its entirety.

### FIELD OF THE INVENTION

10 The present invention relates, in general, to compositions, methods and  
diagnostic kits useful for the detection of fish-killing raphidophyte algae using  
rRNA-targeted probes.

### BACKGROUND OF THE INVENTION

15 Raphidophytes are algae of the class Raphidophyceae. Raphidophytes such  
as *Heterosigma akashiwo* (Hada) ex Sournia and *Fibrocapsa japonica* Toriumi &  
Takano are well known in temperate seas as causative agents for mass finfish kills  
in seapen aquaculture. For example, during the summer of 1989, a large bloom of  
*Heterosigma akashiwo* in Big Glory Bay, Stewart Island, New Zealand, caused  
extensive mortality of caged Quinns salmon (*Onchorynchus tshawytscha*), the loss  
being valued at NZD \$4.5 million [Chang F.H., Anderson, C. & Boustead, N.C.  
20 1990. This event was the first record of a *Heterosigma* (Raphidophyceae) bloom  
with associated mortality of cage-reared salmon in Big Glory Bay, New Zealand.  
*NZ J. Mar. Freshwater. Res.* **24**: 461-469; MacKenzie L. 1991.

Harmful algal bloom research and monitoring has traditionally been based  
on ecological and microbiological measurements which are laborious, time-

consuming, and reliant on experienced operators. The rapid identification and enumeration of harmful raphidophyte species is crucial for the management of cultured finfish, shellfish and wild resources in order to avoid stock loss.

Thus, there is a need to develop a test system for rapid, sensitive and cost effective analysis of Raphidophytes that permits as near as possible real time monitoring of the algae.

### SUMMARY OF THE INVENTION

In order to meet these needs, the present invention is directed to compositions, methods and diagnostic kits useful for the detection of fish-killing raphidophyte algae using rRNA targeted probes. The probes comprise a segment of nucleic acid capable of selectively hybridizing, under selective hybridizing conditions, to large-subunit ribosomal RNA from raphidophytes.

The probes include those oligonucleotide probes having sequences selected from SEQ ID NO: 3 through SEQ ID NO: 23 and homologous sequences. The probes may be utilized in various combinations including pairwise. The probes of the invention may be of the formula  $[X-Y-Z]_n$  where X is a sequence of 0 to 100 nucleotides or nucleotide analogs that are non-homologous to conserved or nonconserved regions of raphidophyte nucleic acid. In the formula, Y is a sequence of 10 to 100 nucleotides or nucleotide analogs that are capable of hybridizing under hybridizing conditions to hypervariable regions of the ribosomal RNA of raphidophytes. Such sequences for Y include those sequences selected from SEQ ID NO:3 through SEQ ID NO:23 and homologous sequences. Furthermore, Z is a sequence of 0 to 100 nucleotides or nucleotide analogs that are non homologous to conserved or non conserved regions of raphidophyte nucleic acid. The sequence of Z may be the same or different from X. Finally, n is 1 to 500 or more.

In addition to compositions, methods are disclosed for the detection of raphidophytes from a marine sample using fluorescent *in situ* hybridization (F.I.S.H.) and sandwich hybridization assays (S.H.A.). These methods comprise the steps of: permeabilizing the species of raphidophyte to be assayed to expose the  
5 ribosomal RNA; contacting the exposed ribosomal RNA, under hybridizing conditions, with oligonucleotide probes capable of selectively hybridizing to the hypervariable regions of the ribosomal RNA of at least one species of raphidophyte; and detecting hybridization complexes as an indication of the presence of the raphidophyte cell in the sample.

10 In addition to compositions and methods, there are disclosed herein diagnostic kits for use in determining the presence of raphidophytes which comprise a synthetic oligonucleotide probe complementary to the aforementioned hypervariable or conserved regions of the ribosomal RNA of a raphidophyte species from a marine sample. The kits may include hybridization buffers.

15 The present invention has utility in providing an easy, sensitive, and specific test for algae which may kill finfish and invertebrates.

### **BRIEF DESCRIPTION OF THE DRAWINGS**

This invention will be better understood by reference to the drawings, in which:

20 Figure 1 shows specificity and isolate comparison for the sandwich hybridization assay probe combinations targeted at *H. akashiwo* and *F. japonica*. In A) the probe combination Het1.25aS/Raphid1F, which is targeted at *H. akashiwo*, was screened against a number of *H. akashiwo* isolates and other raphidophyte species. *H. akashiwo* 1 represents the isolate *H. akashiwo* CAWRO4;  
25 *H. akashiwo* 2 = CAWR05; *H. akashiwo* 3 = CAWR09 and *H. akashiwo* 4 = CAWR14. The *F. japonica* isolate was CAWR02. In B) the probe combination Fib1.25aS/Raphid1F, which is targeted at *F. japonica*, was screened against two

*F. japonica* isolates and other raphidophyte species. *F. japonica* 1 represents the isolate CAWR02 and *F. japonica* 2 represents the isolate CAWR03. The *H. akashiwo* isolate was CAWR05.

5 Figure 2 shows serial dilution standard curves for *F. japonica* and *H. akashiwo* based on the sandwich hybridization assay. In A) a standard curve at 450 nm for *F. japonica* cells spiked into field and filtered sea water background is shown. In B) a standard curve at 450 nm for *H. akashiwo* cells spiked into field and filtered sea water background is shown.

10 Figure 3 shows ribosomal RNA levels as measured using the sandwich hybridization assay over an entire growth cycle for *H. akashiwo* and *F. japonica*. In the graphs B and D the open box data points represent signal at 450 nm, whereas the closed circles represent signal at 655 nm. In A) changes in cell density over a growth cycle for a batch culture of *H. akashiwo* is shown. In B) signal from the sandwich hybridization assay for a constant number of cells over a growth cycle of  
15 *H. akashiwo* is shown. In C) changes in cell density over a growth cycle for a batch culture of *F. japonica* is shown. In D) signal from the sandwich hybridization assay for a constant number of cells over a growth cycle of *F. japonica* is shown.

#### **BRIEF DESCRIPTION OF THE SEQUENCE LISTING**

20 This invention will be better understood by reference to Sequence Listing, in which:

SEQ ID NO: 1 is a Fluorescent *In Situ* Hybridization (F.I.S.H.) positive control probe designated UniC.

SEQ ID NO: 2 is a F.I.S.H. negative control probe designated UniR.

25 SEQ ID NO: 3 is a F.I.S.H. *H. akashiwo* probe designated Het 1.25F.

SEQ ID NO: 4 is a F.I.S.H. *H. akashiwo* probe designated Het 1.5F.

SEQ ID NO: 5 is a F.I.S.H. *H. akashiwo* probe designated Het 2aF.

SEQ ID NO: 6 is a F.I.S.H. *H. akashiwo* probe designated Het 3F.

SEQ ID NO: 7 is a F.I.S.H. *H. akashiwo* probe designated Het. sig2-3'F.

5 SEQ ID NO: 8 is a F.I.S.H. Raphidophyceae probe designated Raphid1F.

SEQ ID NO: 9 is a F.I.S.H. Raphidophyceae probe designated  
Raphid2F.Raphid 2F.

SEQ ID NO. 10 is a F.I.S.H. *F. japonica* probe designated Fib1.25aF.

SEQ ID NO: 11 is a F.I.S.H. *F. japonica* probe designated Fib 1.5F.

10 SEQ ID NO: 12 is a F.I.S.H. *F. japonica* probe designated Fib 2F.

SEQ ID NO: 13 is a F.I.S.H. *F. japonica* probe designated Fib 3F.

SEQ ID NO: 14 is a F.I.S.H. *F. japonica* probe designated Fib.sig2-3'F.

SEQ ID NO: 15 is a Sandwich Hybridization Assay (S.H.A.) *H. akashiwo*  
probe designated Het 1.25 aS.

15 SEQ ID NO: 16 is a S.H.A. *H. akashiwo* probe designated Het 1.25 bS.

SEQ ID NO: 17 is a S.H.A. *H. akashiwo* probe designated Het 3S.

SEQ ID NO: 18 is a S.H.A. *H. akashiwo* probe designated Het.sig2-3'F.

SEQ ID NO: 19 is a Raphidophyceae probe designated Raphid 1F.

SEQ ID NO: 20 is a Raphidophyceae probe designated Raphid 2F.

20 SEQ ID NO: 21 is a *F. japonica* probe designated Fib1.25aS.

SEQ ID NO: 22 is a *F. japonica* probe designated Fib.sig3F.

SEQ ID NO: 23 is a *Chattonella antiqua* probe designated Chat1.25F/S.

SEQ ID NO: 24 is a *Chattonella antiqua* 'D1, D2' domain LSU rDNA gene sequence.

5 SEQ ID NO: 25 is a *Chattonella subsalsa* 'D1, D2' domain LSU rDNA gene sequence.

SEQ ID NO: 26 is a *Fibrocapsa japonica* 'D1, D2' domain LSU rDNA gene sequence.

10 SEQ ID NO: 27 is a *Heterosigma akashiwo* 'D1, D2' domain LSU rDNA gene sequence.

SEQ ID NO: 28 is a *Vacuolaria virescens* 'D1, D2' domain LSU rRNA gene sequence.

SEQ ID NO: 29 is a PCR primer used to amplify the 'D1, D2' domain LSU rRNA gene designated DIR.

15 SEQ ID NO: 30 is a PCR primer used to amplify the 'D1, D2' domain LSU rRNA gene designated D2C.

### DETAILED DESCRIPTION OF THE INVENTION

20 Unless defined otherwise herein, all technical and scientific terms used herein have the same meaning as commonly understood by one of ordinary skill in the art to which this invention belongs. Singleton et al., *Dictionary of Microbiology and Molecular Biology*, second edition, John Wiley and Sons, NY (1994), and Hale and Marham, *The Harper Collins Dictionary of Biology*, Harper Perennial, NY (1991) provide one of skill with a general dictionary of many of the terms used in this invention. Although any methods and materials similar or

equivalent to those described herein can be used in the practice or testing of the present invention, the preferred methods and materials are described. In order to more clearly understand the invention, the following specific definitions are provided:

5           **Raphidophytes:** Raphidophytes are algae of the class Raphidophyceae which includes *Heterosigma akashiwo*, *Fibrocapsa japonica*, *Chattonella antiqua* and others (see Table 1 and C. Van Den Hoek, et. al. Algae, an introduction to phycology, Chapter 10 Heterokontophyta: Class Raphidophyceae, Cambridge University Press, Cambridge, 1995).

10           **Oligonucleotide Probe:** Oligonucleotide probes or polynucleotide probes include both double stranded and single stranded DNA or RNA. The probes may be synthesized synthetically or be recombinantly derived sequences. The probes include a specific oligonucleotide sequence and its complement sequence, e.g. A-T/U and C-G. Oligonucleotide probes include a single oligonucleotide sequence or  
15 a mixture of the given sequences, or a mixture of probes that may contain as a part of the probes single or multiple copies of the given oligonucleotide sequences.

**Probe Compositions:** Probe compositions include probes complementary to raphidophyte rRNA. The probes may be in a pure state or in combination with other probes. In addition, the probes may be in combination with salts or buffers,  
20 and may be in a dried state, in an alcohol solution as a precipitate, or in an aqueous solution. The probes may be a mixture of different probes capable of detecting a single species or two or more species, a mixture of different probes wherein the probes are each able to detect one or more species.

**Open Regions:** Open regions are regions of RNA which have minimal  
25 secondary or tertiary interactions with adjacent nucleotides.

**Closed Regions:** Closed regions are regions of RNA with significant secondary or tertiary interactions with adjacent nucleotides.

**Homologous Sequences:** Homologous sequences are sequences which have sufficient identity to another sequence such that under standard hybridization conditions of moderate stringency the percent hybridization can be shown to exceed 50% of the hybridization between perfectly complementary nucleic acid fragments. Homologous sequences for hybridization complexes between pairs of nucleotides.

**Non Homologous Sequences:** Non homologous sequences are sequences which have sufficient differences to another sequence that under standard hybridization conditions of moderate stringency the percent hybridization is less than 50% of the hybridization between perfectly complementary nucleic acid fragments.

**Marine Sample:** A marine sample is a specimen of sea water or of an organism living within the sea. The term also encompasses a digestive tract specimen from an organism not necessarily living within, but taking nutrition exclusively from, the sea.

**Permeabilize:** Permeabilize means to disrupt a cell to allow for intracellular or extracellular hybridization between nucleic acid probes and rRNA. A permeabilized cell may be lysed to allow extracellular release of rRNA, fixed so as to allow entry of probes into the cell, or both. rRNA is exposed when the rRNA is rendered accessible to hybridization by a complementary segment of nucleic acid.

**Nucleotides:** Nucleotides refer to deoxyribonucleotide or ribonucleotides. These nucleotides may be in the form of a polymer in either single- or double-stranded form as a nucleic acid, and unless otherwise limited, encompass known analogues of natural nucleotides that hybridize to nucleic acids in a manner similar to naturally occurring nucleotides.




**Subsequence:** The subsequence of a particular nucleic acid or polypeptide sequence refers to a region of the nucleic acid or polypeptide smaller than the particular nucleic acid or polypeptide.

**Hybridization Wash Conditions:** Hybridization wash conditions refer to those wash conditions for nucleic acid hybridization experiments such as Southern and northern hybridizations which are sequence dependent. In nucleic hybridization experiments, hybridization complexes are formed between homologous sequences. An extensive guide to the hybridization of nucleic acids is found in Tijssen, *Laboratory Techniques in Biochemistry and Molecular Biology -- Hybridization with Nucleic Acid Probes*, Part I, Chapter 2 "Overview of principles of hybridization and the strategy of nucleic acid probe assays", Elsevier, New York (1993) which is hereby incorporated by reference. Stringency of the wash conditions are dependent on numerous factors including the temperature, components and other factors. Generally, highly stringent wash conditions are selected to be about 5°C lower than the thermal melting point ( $T_m$ ) for the specific sequence at a defined ionic strength and pH. The  $T_m$  is the temperature (under defined ionic strength and pH) at which 50% of the target sequence hybridizes to a perfectly matched probe. Very stringent conditions are selected to be equal to the  $T_m$  point for a particular probe.

Taking into account the above definitions, the present invention relates to compositions, methods and diagnostic kits for the detection of fish killing raphidophyte algae using RNA targeted probes. The probes comprise a segment of nucleic acid capable of selectively hybridizing, under selective hybridization conditions, to large-subunit ribosomal RNA from raphidophytes. The probes include those sequences selected from SEQ ID NO: 3 through SEQ ID NO: 23 and homologous sequences.

1. **Probes**

The probes of the present invention are generally of the formula:

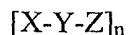

In the formula, X is a sequence of 0 to 100 nucleotides or nucleotide  
5 analogs that are non-homologous to conserved or non-conserved regions of  
raphidophyte nucleic acid. Furthermore, Y is a sequence of 10 to 100 nucleotides  
or nucleotide analogs that are capable of hybridizing under hybridizing conditions  
to hypervariable regions of the ribosomal RNA of raphidophytes. Such nucleotides  
or nucleotide analogs include those sequences selected from SEQ ID NO: 3  
10 through SEQ ID NO: 23.

In the formula, Z is a sequence of 0 to 100 nucleotides or nucleotide analogs  
that are non homologous to conserved or non conserved regions of raphidophyte  
nucleic acid. The sequence of Z may be the same or different from X. Finally, n is  
1-500 or more. In the probes where n is greater than 1, Y can be the same or  
15 different sequences of nucleotides having hybridization capability. The probe can  
be free or contained within a vector sequence (e.g., plasmids, viruses or cosmid).

The nucleic acid sequence of the claimed probes include homologous  
synthetically derived or recombinant nucleic acid sequences which have sufficient  
identity with the claimed sequences that they substantially hybridize with regions  
20 complementary to the claimed probes to form hybridization complexes. By  
"substantially" it is meant that under standard hybridization conditions of moderate  
stringency, percent hybridization can be shown to exceed 50% of the hybridization  
between perfectly complementary nucleic acid fragments.

The probes of the present invention substantially bind under selective  
25 hybridizing conditions to regions of raphidophyte rRNA having minimal secondary  
or tertiary interactions with adjacent nucleotides known as open regions. By  
"substantially bind" it is meant that the probes do not comprise significant

sequences that bind to regions that are available for hybridization only after heating, that is, regions with significant secondary and tertiary structure (closed regions). In practical terms, such probes will generally not comprise any more than 10 flanking nucleotides (either 5' or 3') which would bind to closed regions. More specifically, compositions of polynucleotide probes complementary to open regions are claimed that are complementary to either the hypervariable or, alternatively, to the conserved regions of rRNA of raphidophyte species.

The probes of the present invention may be compounds of RNA, DNA, or RNA/DNA chimeras. In the probes, analogs of nucleotides may be substituted for naturally occurring nucleosides. Nucleotide analogs offer greater stability, resistance to nuclease activity and ease of signal attachment. The term "probe" is intended to embrace all functionally equivalent sequences. Equivalent raphidophyte probes may also consist of the given sequence, concatemers of the sequence, or probes flanked by about 10 or less bases of any degree of complementarity to the native sequences flanking the complementary region of raphidophyte rRNA.

## 2. Probes Complementary to rRNA of Raphidophytes

The degree of complementarity (homology) required for detectable binding of probes of the present invention with the rRNA of raphidophytes will vary in accordance with the stringency of the hybridization medium and/or wash medium. The degree of complementarity will optimally be 100 percent; however, it should be understood that minor variations between the rRNA and the disclosed probes may still provide for selective hybridization to a particular rRNA without undesired cross-hybridization to other accessible nucleic acids in the sample. Such minor variations may be compensated for by reducing the stringency of the hybridization and/or wash medium as described below. Thus, despite the lack of look complementarity under reduced conditions of stringency, functional probes having minor base differences from their rRNA targets are possible. Therefore, under

hybridization conditions of reduced stringency, it may be possible to slightly modify the disclosed probes while maintaining an acceptable degree of specificity to detect the desired raphidophyte species present in the sample.

### 3. Synthesis or Isolation of the Probes

5           The probes of the present invention may be chemically synthesized using commercially available methods and equipment. Methods of synthesizing nucleic acids are well known in the art. Deoxynucleotides may be synthesized chemically according to the solid phase phosphoramidite triester method described by Beaucage and Caruthers, *Tetrahedron Letts.* 22(20):1859-1862 (1981), e.g., using  
10           an automated synthesizer, e.g., as described in Needham-VanDevanter et al., *Nucleic Acids Res.* 12:6159-6168 (1984).

          To obtain large quantities of oligonucleotide probes, one can also clone the desired sequence using traditional cloning methods, such as described in Maniatis, T., et al. *Molecular Cloning: A Laboratory Manual*, Cold Spring Harbor, New  
15           York (1982), or one can produce the probes by chemical synthesis using commercially available DNA synthesizers. An example of cloning would involve insertion of the cDNA for the ribosomal RNA into a replication vector, such as pBR322, M13, or into a vector containing the SP6 promotor (e.g., generation of  
20           single-stranded RNA using SP6 RNA polymerase), and transformation of a bacterial host. The DNA probes can be purified from the host cell by lysis and nucleic acid extraction, treatment with selected restriction enzymes, and further isolation by gel electrophoresis. The use of polymerase chain reaction (PCR) technology can also be used to obtain large quantities of probe. (See U.S. Patent No. 4,683,202.)

### 25           4. Uses of the probes

          The probes of the invention are useful for detecting raphidophytes in a marine sample. The present invention is thus also directed to methods of detecting

a raphidophyte species from a marine sample. These methods include the steps of:  
(a) permeabilizing the cells of the raphidophyte species to expose the ribosomal  
RNA; (b) contacting the ribosomal RNA under hybridizing conditions with  
polynucleotide probes capable of selectively hybridizing to a hypervariable or  
5 conserved region of the ribosomal RNA of species; and (c) detecting hybridization  
complexes as an indication of the presence of the species in the sample.

Marine samples for use in this invention can be obtained by any number of  
methods well known to the skilled artisan including tow samples of marine waters.  
The samples may be collected using an aquatic autosampler as described in  
10 applicant's co-pending United States patent application Serial Number 09/319,333  
entitled "Aquatic Autosampler" which is hereby incorporated by reference. The  
samples may be subsequently processed, for example, to remove precipitated  
material, or gently filtered to concentrate the sample or exclude organisms of  
particular size, or cultured to enrich or deplete the population of particular  
15 organisms. Conveniently, the sample is dispersed in a buffer protective of, or  
compatible with rRNA [Boney, A.D., *Phytoplankton, Studies in Biology*, no. 52  
(1979)].

The sample, or a portion thereof, may be permeabilized for use in a  
sandwich hybridization assay in a lysis buffer such as disclosed in Van Ness et al.,  
20 *Nucl. Acids. Res.* **19**:5143-5151 (1991), and PCT application WO 93/24659, both  
incorporated herein by reference. Lysing solutions are well known in the art and  
are typically composed of a buffered detergent solution having a divalent metal  
chelator or a buffered chaotropic salt solution containing a detergent (such as  
SDS), a reducing agent and a divalent metal chelator (EDTA). Generally, these  
25 buffers are between pH 7.0 and 9.0, and contain both chelating agents and  
surfactants.

Mechanical methods, including French press, nitrogen cavitation, bead  
beater, ultrasound sonification, and heating, may also be employed to permeabilize

the cell. Alternatively, samples may be collected and dispersed in a lysing solution that also functions as a hybridization solution, such as 3M guanidinium thiocyanate (GuSCN), 50mM Tris (pH 7-6), 10mM EDTA, 0.1% sodium dodecylsulfate (SDS), and 1% mercaptoethanol [Maniatis, *et al.*, *Molecular Cloning: A Laboratory*  
5 *Manual*, Cold Spring Harbor, NY (1982)].

### 5. Hybridization Assays

The probes of the invention may be utilized to assay for raphidophytes using various hybridization assays. Assay test protocols for use in this invention are those of convention in the field of nucleic acid hybridization. Assay test  
10 protocols include both single phase hybridizations, where the target and probe polynucleic acids are both in solution and mixed phase hybridizations, where either the target or probe polynucleotides are fixed to an immobile support. Mixed phase hybridizations include non-sandwich type assays. Whole cell hybridizations may also be employed using methods well known in the art and exemplified herein. The  
15 assay test protocols known to the skilled artisan are varied and are not to be considered a limitation of this invention.

Various hybridization solutions may be employed. Hybridization solutions generically include from about 20 to 60% volume, preferably 30%, of a polar organic solvent. A common hybridization solution employs about 5% v/v  
20 formamide, about 0.5 to 1M sodium chloride, about 0.05 to 0.1M buffers, such as sodium citrate, Tris-HCl, PIPES or HEPES (pH range about 6-9), about 0.05 to 0.2% detergent, such as sodium dodecylsulfate, or between 0.5-20mM EDTA, 0.01-0.051% Ficoll (about 300-500 kilodaltons), 0.01-0.05% polyvinylpyrrolidone (about 250-500 KDa), and 0.01-0.05% serum albumin. Also included in the typical  
25 hybridization solution will be unlabeled carrier nucleic acids from about 0.1 to 5 mg/ml, fragmented nucleic DNA, e.g., calf thymus or salmon sperm DNA, or yeast RNA, and optionally from about 0.5 to 2% wt./vol. glycine. Other additives may also be included, such as volume exclusion agents which include a variety of polar

water-soluble or swellable agents, such as polyethylene glycol, anionic polymers such as polyacrylate or polymethylacrylate, or polystyrene sulfonic acid and anionic saccharidic polymers, such as dextran sulfate.

5 An alternative hybridization solution may be employed including about 2 to 4M GuSCN, preferably 3M, about 0.01 to 0.1M Tris (pH range about 6.0 to 8.9), a detergent such as sodium dodecyl sulfate in concentrations of about 0.1 to 5% (wt./vol.), and about 0.01 to 0.1M EDTA. Other additives may also be included such as carrier DNA or RNA, or protein such as bovine serum albumin or gelatin. Stringency of the hybridization solution can be adjusted by the addition of about 0  
10 to 10% formamide, usually 5%.

The particular hybridization technique utilized is not essential to the invention. Hybridization techniques are generally described in *Nucleic Acid Hybridization: A Practical Approach*, Ed. Hames, B.D. and Higgins, S.J., IRL Press (1987); Gall *et al.*, *Proc. Natl. Acad. Sci., U.S.A.*, **63**:378-383 (1969), and  
15 John *et al.* *Nature*, **223**:582-587 (1969). As improvements are made in hybridization techniques, they can readily be applied.

Regardless of the assay test protocol being used, the raphidophyte cells or cell contents are to remain in contact with a hybridization solution for an extended period of time. In single phase assays, the double-stranded duplexes may be  
20 separated from single-stranded nucleic acid by S<sub>1</sub> nuclease digestion followed by precipitation of duplex molecules, or by selective binding to hydroxyapatite. In mixed phase assays, the support-immobilized nucleic acid is typically introduced into a wash solution having analogous concentrations of sodium chloride, buffers, and detergent, as provided in the hybridization solution. The time period for which  
25 the support is maintained in the wash solution may vary from several minutes to three hours or more.

The amount of labeled probe which is added to the hybridization solution may vary widely, depending upon the nature of the label, the amount of the labeled

probe which can reasonably bind to the cellular target nucleic acid, and the stringency of the hybridization medium and/or wash medium. Generally, substantial excesses of probe over the stoichiometric amount of the target nucleic acid will be employed to enhance the rate of binding of the probe to the target DNA.

Various degrees of stringency of hybridization can be employed. As the conditions for hybridization become more stringent, there must be a greater degree of complementarity between the probe and the target for duplex formation to occur. The degree of stringency can be controlled by temperature, ionic strength, pH and the presence of a partially denaturing solvent such as formamide. For example, the stringency of hybridization is conveniently varied by changing the polarity of the reactant solution through manipulation of the concentration of formamide within the range of 0% to 50%.

Either the hybridization or the wash medium can be stringent. Typically, for mixed phase assays, it is the wash solution that most often determines the stringency and facilitates dissociation of mismatched duplexes. After rinsing the support at room temperature with a dilute buffered sodium chloride solution, the support may now be assayed for the presence of duplexes in accordance with the nature of the label.

In F.I.S.H. assays cells are fixed, treated with hybridization buffer and then hybridized with probe as described in the Example section below. The cells are then analyzed by epifluorescence microscopy.

In a sandwich-type assay a primary component is a solid support. The solid support has absorbed to it or covalently coupled to it immobilized nucleic acid probe that is unlabeled and complementary to one portion of the rRNA sequence. Probes hybridize to regions of the ribosomal RNA with minimal secondary and tertiary interactions. The advantage of such probes is that the hybridization can be carried out without the additional step of heat denaturing the sample nucleic acid.



For the sandwich-type assay, the test sample suspected of containing one or more raphidophytes is contacted with the solid support in a hybridization medium. In the procedure, a second soluble-labeled probe complementary to a different sequence of the rRNA of the raphidophyte is hybridized to the rRNA that has formed a hybridization duplex with the immobilized nucleic acid probe on the solid support. A probe to a hypervariable region and a probe to a conserved region of rRNA of raphidophyte may each function as either a capture or signal probe. The entire assay takes place at 30°C.

The presence of raphidophyte species assayed from the marine sample is then determined in accordance with the label being used. It should be noted that in sandwich hybridization assays the second probe can be added simultaneously with the test sample to the hybridization assay. In addition, the second probe can hybridize to either a conserved or to a hypervariable region of the rRNA.

Compared to the whole cell F.I.S.H. method, sandwich hybridization is several fold faster and far less taxing on the individual executing the protocol [Scholin, C. A., et al. DNA probes and a receptor-binding assay for detection of *Pseudo-nitzschia* (Bacillariophyceae) species and domoic acid activity in cultured and natural samples. J. Phycol. 35: 1356-1367 (1999)]. Equipment used to apply this method includes a standard filtration manifold, heating block, micro pipette, and a robotic processor. The latter costs several thousand dollars less than an average epifluorescence microscope and is available commercially. It is possible to apply this method outside of a laboratory. Using current technology, the lower limit of detection (LLD) for sandwich hybridization is roughly  $2.5 - 5 \times 10^2$  raphidophyte cells  $0.2 \text{ ml}^{-1}$  of lysate. For whole water samples, it is possible to detect as few as several hundred to  $10^3$  cells  $\text{l}^{-1}$  using the current sandwich hybridization assay and instrumentation described above.

6. Detection of Raphidophytes

Various detection labels may be utilized. Where the label is radioactive, the presence of probe can be detected in a scintillation counter. More conveniently, in mixed phase assays, the substrate can be dried and exposed to X-ray film in any number of conventional autoradiographic protocols. Autoradiographic detection is typically employed with  $^3\text{H}$ ,  $^{125}\text{I}$ ,  $^{35}\text{S}$ ,  $^{14}\text{C}$ , or  $^{32}\text{P}$  labeled probes or the like. The choice of radioactive isotope depends on research preferences due to ease of synthesis, stability and half lives of the selected isotopes.

Where the label is fluorescent, the sample is detected by first irradiating it with light of a particular wavelength. The sample absorbs this light and then emits light of a different wavelength which is picked up by a detector.

Where the label is a hapten or antigen, the sample can be detected by using antibodies. In these systems, a signal is generated by attaching fluorescent or enzyme molecules to the antibodies; in some cases the antibody is labeled with a radioactive probe as described in Tijssen, P., *Practice and Theory of Enzyme Immunoassays, Laboratory Techniques in Biochemistry and Molecular Biology*, Burdon, R.H., van Knippenberg, Ph.H., Eds., Elsevier, pp. 9-20 (1985), which is hereby incorporated by reference.

One method of detection is enzymatic detection in conjunction with biotin. Although fluorescence is an alternative label, enzymatic labels, in combination with avidin or streptavidin such as biotinylated peroxidase or alkaline phosphatase, are preferred. Enzyme-conjugated avidin or streptavidin can also be used to directly bind the enzyme to the probe. Preferred enzymes are peroxidase or alkaline phosphatase.

Non-radioactive probes are often labeled by indirect means. For example, a ligand molecule is covalently bound to the probe. The ligand then binds to an anti-ligand molecule which is either inherently detectable or covalently bound to a

detectable signal system, such as an enzyme, a fluorophore, or a chemiluminescent compound. Ligands and anti-ligands may be varied widely. Where a ligand has a natural anti-ligand, namely ligands such as biotin, thyroxine, and cortisol, it can be used in conjunction with its labeled, naturally occurring anti-ligands. Alternatively, any haptenic or antigenic compound can be used in combination with an antibody.

Probes can also be labeled by direct conjugation with a label. For example, cloned DNA probes have been coupled directly to horseradish peroxidase or alkaline phosphatase, [Renz, *et al.*, *Nuc. Acids Res.* **12**:3435-3444 (1984)] and synthetic oligonucleotides have been coupled directly with alkaline phosphatase [Jablonski *et al.*, *Nuc. Acids. Res.* **14**:6115-6128 (1986)]. A general reference for various detection methods can be found in Hames, B. D. and Higgins, S. J., *Nucleic Acid Hybridization*, IRL Press, Oxford (1985). References for sandwich assay with DNA probes are Dunn and Hassell, *Cell*, Vol. **12**, pp. 23-26 (1977), and Ranki, *et al.*, U.S. Pat. No. 4,486,539. All of these articles are hereby incorporated by reference.

## 7. Kits

The oligonucleotide or polynucleotide acid probes of this invention can be included in a kit which can be used to rapidly determine the presence or absence of raphidophyte species from a marine sample. The kit includes all components necessary to assay for the presence of these species. The kit includes a stable preparation of rRNA probes, hybridization solution in either dry or liquid form for the hybridization of target and probe polynucleotides, as well as a solution for washing and removing undesirable and nonduplexed polynucleotides, and a substrate for detecting the labeled duplex. The probes may be labeled or unlabeled. The kit will also include directions for using the probes.

A more specific embodiment of this invention embraces a kit that utilizes the concept of the sandwich hybridization assay (S.H.A.). This kit would include as a first component, vials for containment of a marine sample and buffers for the

permeabilization of the sample. A second component would include media in either dry or liquid form for the hybridization of target and probe polynucleotides, as well as for the removal of undesirable and nonduplexed forms by washing. A third component includes a solid support (such as a dipstick) upon which is fixed or to which is conjugated unlabeled nucleic acid probe(s) that is (are) complementary to a part of the rRNA of the species assayed for. In the case of multiple target analysis more than one capture probe, each specific for its own ribosomal RNA, will be applied to different discrete regions of the dipstick. A fourth component would contain a labeled or unlabeled probe that is complementary to a second and different region (conserved or hypervariable) of the same rRNA strand to which the immobilized, unlabeled nucleic acid probe of the third component is hybridized.

The probe components described herein include coordinations of probes in dry form, such as lyophilized nucleic acid or in precipitated form, such as alcohol precipitated nucleic acid or in buffered solutions. The various reagents for the detection of labeled probes and other miscellaneous materials for the kit, such as instructions, positive and negative controls, and containers for conducting, mixing, and reacting the various components, would complete the assay kit. Such kits would include instruction cards and vials containing the various solutions necessary to conduct a nucleic acid hybridization assay. These solutions would include lysing solutions, hybridization solutions, combination lysing and hybridization solutions, and wash solutions. The kits would also include labelled probes. Standard references for comparison of results may also be provided for an easy estimate of the number of raphidophyte species in a given solution. Depending upon the label used additional components may be needed for the kit, e.g., enzyme labels require substrates.

From the foregoing description it will be clear to those of skill in the art that the disclosed probes are useful for identifying raphidophytes collected from pure cultures or nature. Moreover, the technique offers promising means by which one could quantify these species rapidly. The protocol is neither extremely complicated

or demanding, and with minimal training all individuals should be able to execute both methods.

The invention is now illustrated by way of the following non-limiting examples.

5

### Examples

#### **Example 1**

##### Algae Cultures

10 The algae cultures utilized to isolate the probes of this invention are listed in Table 1. All cultures are currently held at the Cawthron Institute, Nelson, New Zealand (K. Ponikla). Cultures are maintained under a 12:12 light/dark photoperiod at 20°C in f/2 media (Guillard, R. R. L. 1975. Culture of phytoplankton for feeding marine invertebrates. In: Culture Of Marine Invertebrate Animals (ed. by W. L. Smith & M. H. Chanley), pp. 29-60, Plenum Press, New York.).

**Table 1. Species probed in this study and culture collection numbers.**

Species	Class	Culture collection number
<i>Chattonella antiqua</i>	Raphidophyceae	CAWRO1
<i>Chattonella subsalsa</i>	Raphidophyceae	CCMP217
<i>Fibrocapsa japonica</i>	Raphidophyceae	CAWRO2, O3
<i>Heterosigma akashiwo</i>	Raphidophyceae	CAWRO4, O5, O6, O9, 14(ABC2)
<i>Heterosigma</i> sp.	Raphidophyceae	CAWR10
<i>Olisthodiscus luteus</i>	Uncertae sedis	NIES-15
<i>Nannochloropsis oculata</i>	Eustigmatophyceae	LB2164

\*Notes on Table 1:

- CAW -Cawthron Institute, Nelson, New Zealand, Algal Collection  
-CAW14(ABC<sup>2</sup>) was provided by Dr. Rita Horner,  
University of Washington, Seattle, Washington, USA.
- CCMP Provasoli-Guillard National Centre for Culture of Marine  
Phytoplankton, West Boothbay Harbor, Maine, USA
- LB -The University of Texas of Austin, Austin, Texas, USA, Culture Collection of Algae
- NIES -National Institute for Environmental Studies, Japan, Microbial Culture  
Collection, City, Japan

## Example 2

### DNA Extraction

Cells from Example 1 grown to mid-log phase growth were harvested by centrifugation. Approximately 10mg of algal cells was added to 700 $\mu$ L of CTAB buffer (50mM Tris pH 8.0, 0.7M NaCl, 10mM EDTA, 1% CTAB). 1 $\mu$ L of 0.1% of  $\beta$ -mercapto-ethanol was then added and the mixture vortexed to resuspend the pellet. After the pellet had been resuspended 30 $\mu$ L of 20mg.mL<sup>-1</sup> Proteinase K was added, gently mixed, then incubated at 60°C for 1 hour. If the pellet was still intact after 1 hour a further 5 $\mu$ L of Proteinase K was added and incubated for another hour. After the incubation an equal volume of 24:1 chloroform:isoamyl alcohol was added, vortexed and placed on a shaker for 5 minutes. The mixture was then centrifuged at 13,000 rpm for 5 minutes at 4°C and the supernatant transferred to a new tube. An equal volume of 25:24:1 phenol:chloroform:isoamyl alcohol was added, vortexed and placed on a shaker for 5 minutes. The mixture was then centrifuged at 13,000 rpm for 5 minutes at 4°C and the supernatant transferred to a new tube. This step was continued until the interphase was clear. Another 24:1 chloroform:isoamyl alcohol extraction was then done as for the first step and the supernatant transferred to a new tube. An equal volume of isopropanol was added to the supernatant, mixed and placed at -20°C for at least an hour. The DNA was pelleted by centrifugation at 13,000 rpm at 4°C for 30 minutes, then washed 2 times with 70% ethanol, using a 5 minute centrifugation step at 13,000 rpm at 4°C. The DNA pellet was dried under vacuum and resuspended in 50 $\mu$ L of ddH<sub>2</sub>O.

### Example 3

#### PCR Amplification

Amplification of DNA from Example 2 was achieved using the polymerase chain reaction (PCR). The following primers were used to amplify the 'D1, D2' domain LSU rDNA gene:

D1R (forward) 5' ACC CGC TGA ATT TAA GCA TA 3' (SEQ ID NO: 29)

D2C (reverse) 5' CCT TGG TCC GTG TTT CAA GA 3' (SEQ ID NO: 30)

The primers are targeted at conserved regions at positions 24-45 for D1R and 733-714 for D2C, relative to the *Prorocentrum micans* LSU rRNA gene sequence.

Approximately 10ng of genomic DNA was added to a PCR mix containing 0.25mM dATP, 0.25mM dCTP, 0.25mM dGTP, 0.25mM dTTP, 2.5mM MgCl<sub>2</sub>, 50mM KCl, 50ng of each primer, 1.25 units of *Taq* polymerase (AmpliTaq™) polymerase and sterile water to a final volume of 50μL. A DNA thermal cycler was used to subject the reaction to an initial cycle of denaturation (Hot Start, 94°C for 3 minute), and then 15 cycles of denaturation (94°C for 30 seconds), annealing (60°C each round for 30 seconds) and extension (72°C for 1 minute). These initial 15 cycles were followed by 15 more rounds of denaturation (94°C for 30 seconds), annealing (55°C for 30 seconds) and extension (72°C for 1 minutes) and a final extension of 7 minutes. To check PCR fragment size and purity, 5μL of PCR product was run out on a 0.8% agarose gel, stained with ethidium bromide and viewed on a UV transilluminator.

### Example 4

#### Automated DNA sequencing

Sequencing reactions were carried out using an ABI 373A 'stretch' automated sequencer. Both the coding and non-coding strands were sequenced from a pooled PCR product (at least 2 PCR reactions) from Example 3, using the

PCR primers in dye terminator reactions. The sequences were checked and corrected with the SEQED and gelassemble programs, using the sequence electrophoretogram for visual assurance (Genetics Computer Group. 1994. *Program manual for GCG package, Version 8*. Wisconsin, USA).

5      **Example 5**

**Probe Design and Synthesis**

10      The D1 and D2 domains of the LSU rRNA gene (D1, D2 LSU rRNA) were sequenced (SEQ ID N0:24-28) to provide information for the design of species-specific rRNA-targeted oligonucleotide probes and to construct a phylogenetic framework for intra-class organization of the Raphidophyceae. The rRNA gene transcript contains a mosaic array of conserved and hypervariable domains, which provide unique nucleotide 'signatures' that allow the design of oligonucleotide probes which can be diagnostic from Kingdom to species or even strain-specific levels. Oligonucleotide probes generally are targeted at rRNA because the high  
15      copy number present in cells offers a naturally amplified target for probe localization.

20      Species-specific oligonucleotide probes were designed by looking for unique 'signatures' in the alignment and aiming for a 50% GC content with a probe/target melting temperature ( $T_m$ ) of between 60°C and 75°C depending on the application. The initial oligonucleotide probes were then analysed by the program OligoTech, version 1.0 (Oligos Etc., Eugene, Oregon) to check for intramolecular folding and/or homodimer problems and adjusted as required. These oligonucleotide probes were synthesized commercially with either a fluorescein or biotin moiety coupled to the 5' end (Oligos Etc., Eugene, Oregon, USA). The  
25      probes were resuspended in ddH<sub>2</sub>O at a concentration of 250 ng  $\mu\text{L}^{-1}$ , pipetted to multiple 0.6 mL tubes, then vacuum desiccated and stored at -80°C. Working stocks were resuspended in TE buffer (10 mM Tris HCl, 1 mM EDTA, pH 7.4) at a



final concentration of 250 ng  $\mu\text{L}^{-1}$  for F.I.S.H. probes and 100 ng  $\mu\text{L}^{-1}$  for the S.H.A.

5 An oligonucleotide probe (Chat1.25F/S) has been successfully developed for *Chattonella* spp. which displays good sensitivity and specificity in fluorescent *in situ* hybridization assays. The design of the probe was based on the sequences from the D1, D2 domains of the large subunit rRNA (LSU rRNA) gene for *Chattonella subsalsa* and *C. antiqua* (the *C. antiqua* gene sequence is homologous for the species *C. marina* and *C. ovata*) and an alignment of raphidophyte sequences. The probe Chat1.25F/S is in the same general position within the LSU  
10 rRNA gene as the probes for *Heterosigma akashiwo* (Het1.25aS) and *Fibrocapsa japonica* (Fib1.25S). The addition of an assay for *Chattonella* species completes the coverage of probes for species from the Raphidophyceae which are known to cause mortality of caged fish.

15 The sequences were aligned, revealing several unique species-specific nucleotide sequences and from these unique regions we designed a suite of oligonucleotide probes for *Heterosigma akashiwo*, *Fibrocapsa japonica*, *Chattonella* and raphidophytes in general. The oligonucleotide sequences are shown (SEQ ID NO:1-23) in Table 2.

**Table 2. Oligonucleotide probe sequences, targets and labels**

Probe	Target Species	Sequence and label type (5' - 3')	SEQ ID NO:
UniC	positive control	Fluorescein-GWATTACCGCGGCKGCTG	1
uniR	negative control	Fluorescein-CAGCMGCCGCGGTAAATWC	2
Het1.25F	<i>H. akashiwo</i>	Fluorescein-CGACTGAGCACGCACCTTT	3
Het1.5F	<i>H. akashiwo</i>	Fluorescein-GCGACGGCAAAAAGACCAGGA	4
Het2aF	<i>H. akashiwo</i>	Fluorescein-GCATGTTGAAACGCTCCAG	5
Het3F	<i>H. akashiwo</i>	Fluorescein-AGCAAAGGTCCTCCGTCCTA	6
Het.sig2-3'F	<i>H. akashiwo</i>	Fluorescein- TACTCTCTTTTCAAAGTCTTTTCATC	7
Raphid1F	Raphidophyceae	Fluorescein-CCGCTTCACTCGCCGTTACTAG	8
Raphid2F	Raphidophyceae	Fluorescein- TCATCTTTCCCTCACGGTACTTGTT	9
Fib1.25aF	<i>F. japonica</i>	Fluorescein-CGGCTGGACACGCTTCTGT	10
Fib1.5F	<i>F. japonica</i>	Fluorescein-CAGCACGAAATATGACCCCCG	11
Fib2F	<i>F. japonica</i>	Fluorescein- CCATGGGACACAGCGCGCACTAC	12
Fib3F	<i>F. japonica</i>	Fluorescein-TACAAACCAAGGTGCACTAATG	13
Fib.sig2-3'F	<i>F. japonica</i>	Fluorescein- AACTCTCTTTCCAAAGTTCTTTTCATC	14
Het1.25aS	<i>H. akashiwo</i>	Biotin-ACCACGACTGAGCACGCACCTTT	15
Het1.25bS	<i>H. akashiwo</i>	Biotin-AGCCCGGGACCACGACTGAG	16
Het3S	<i>H. akashiwo</i>	Biotin-GAGCAAAGGTCCTCCGTCCTAAC	17
Het.sig2-3'F	<i>H. akashiwo</i>	Fluorescein- TACTCTCTTTTCAAAGTCTTTTCATC	18
Raphid1F	Raphidophyceae	Fluorescein-CCGCTTCACTCGCCGTTACTAG	19
Raphid2F	Raphidophyceae	Fluorescein- TCATCTTTCCCTCACGGTACTTGTT	20
Fib1.25aS	<i>F. japonica</i>	Biotin-CGGCTGGACACGCTTCTGTAG	21
Fib.sig3F	<i>F. japonica</i>	Fluorescein- AACTCTCTTTCCAAAGTTCTTTTCATC	22
Chat1.25F/S	<i>Chattonella</i> spp.	(Fluorescein/Biotin*)- AGAGTAGCTGAGCACGCATCTCT	23

\*Fluorescein or Biotin label

## Example 6

### Screening of Candidate Probes

Initial screening of probes was conducted at approximately 10°C below the T<sub>m</sub> (melting temperature) and varied according to the degree of specificity required using Fluorescent In Situ Hybridization (F.I.S.H.). This initial hybridization temperature gives a good indication of probe accessibility to target rRNA and the signal intensity. In this study, the specificity of the F.I.S.H. reaction was controlled by the hybridization temperature, not the wash temperature and/or salt concentration.

The relative signal intensity of the species-specific and raphidophyte-specific probes spanning the D1, D2 rRNA gene were compared against a positive control probe (UniC) and a negative control probe (UniR) (SEQ ID No. 1-2). The positive control probe is targeted at universally conserved sequence of the small subunit rRNA gene (SSU rRNA) and should bind to the SSU rRNA of all organisms. The negative control probe is the complement of UniC and does not bind to any known rRNA target and therefore helps to determine the level to which probe retention is non-specific.

Approximately 5-10 mL of mid-exponential culture, net tow or whole (unconcentrated) seawater sample was pipetted gently into a 50 mL polypropylene, conical bottom, disposable centrifuge tube containing an ethanol/saline fixative ((2 mL ddH<sub>2</sub>O, 3 mL 25 x SET buffer (3.75 M NaCl, 25 mM EDTA, 0.5 M Tris HCl, pH 7.8) and 25 mL of 95% ethanol)). The mixture was left to stand at room temperature for at least 15 minutes before gently mixing by inversion, then allowed to stand an additional hour. Aliquots of the samples were filtered onto either 13 mm diameter, 1.2 µm pore size Isopore (Millipore) or Cyclopore (Whatman) polycarbonate membranes in a custom filter tube (Miller P.E. & Scholin C.A. 1998. Identification and enumeration of cultured and wild *Pseudo-nitzschia*

(Bacillariophyceae) using species-specific LSU rRNA-targeted fluorescent probes and filter-based whole cell hybridization. *J. Phycol.* **34**: 371-382). The filtered samples were washed with 1 mL of 5 x SET hybridization buffer (5X SET, 0.1% v/v Igepal, 10  $\mu\text{g}\cdot\text{mL}^{-1}$  poly A, Sigma), then resuspended in 0.5 mL 5 x SET hybridization buffer, to which probe was added at a final concentration of 5  $\text{ng}\cdot\mu\text{L}^{-1}$ . The filter manifold was then transferred to either a dark dry incubator or water bath for 30 min at 50°C. After the incubation, the excess probe was removed by washing with 1 mL of 5 x SET buffer at room temperature for 2 min to remove excess unbound probe. The filters were then removed from the manifold and mounted on glass slides, sample side up, covered with 20  $\mu\text{L}$  of SlowFade Light (Molecular Probes, Eugene, OR) and then mounted with cover slips. Samples were examined by epifluorescence microscopy using a Zeiss Axioskop fitted with a fluorescein band-pass filter set (excitation 465-495 nm; dichoric 505 nm; emission 515-555 nm) and Olympus 10AD 35 mm camera system.

The relative signal intensity of the putative species and raphidophyte-specific LSU rRNA-targeted probes were compared against positive (UniC) and negative (UniR) control probes (Table 3). Cell staining intensity was scored visually using the criteria of Miller and Scholin (Miller P.E. & Scholin C.A. 1996. Identification of cultured *Pseudo-nitzschia* (Bacillariophyceae) using species-specific LSU rRNA-targeted fluorescent probes. *J. of Phycol.* **32**: 646-655.): cells with signal intensity similar to the positive control were scored as ++, signal intensity equivalent to the negative control was scored as --, and signal intensities clearly above the negative but below the positive control as +-. All F.I.S.H. probes were applied to log phase cultures of *H. akashiwo*, *F. japonica* and *Chatonella* spp.

A preliminary screening of oligonucleotide probes using the F.I.S.H. assay was performed and the results are presented in Table 3. The positive control probe

gave a clear and strong signal for all species examined. Furthermore, there was no non-specific retention of the negative control probe.

Of the probes tested, three putative species-specific probes for *Heterosigma akashiwo* (Het1.25F, Het3F and Het. sig.2-3'F) and *Fibrocapsa japonica* (Fib1.25aF, Fib2F and Fib. sig.2-3'F) gave signal which was equivalent to the positive control. The Chat1.25F also labeled its target species with intensity similar to positive control. Both raphidiophyte signal probes (Raphid1F & 2F) gave signals which were comparable to the positive control. From these results, the probes which generated good signal were further tested for specificity against a range of isolates for the target species as well as non-target species. The results for this analysis are presented in Table 4.



Table 4. Specificity and sensitivity of probes using the F.I.S.H. assays.

Species	Probes										
	Positive control	Negative control	Het1.25F	Het3F	Het.sig2-3'F	Fib1.25aF	Fib2F	Fib.sig2-3'F	Chat 1.25F	Raphid1F	Raphid2F
<i>H.akashiwo</i>	++	--	++	++	++	+-	--	+-	--	++	++
<i>Heterosigma</i> sp.	++	--	++	++	++	+-	--	+-	--	++	++
<i>F.japonica</i>	++	--	+-	--	+-	++	++	++	--	++	++
<i>C.antiqua</i>	++	--	+-	--	+-	+-	--	+-	++	++	++
<i>C.subsalsa</i>	++	--	+-	--	+-	+-	--	+-	++	++	++
<i>O.luteus</i>	++	--	--	--	--	--	--	--	--	+-	+-
<i>N.oculata</i>	++	--	--	--	--	--	--	--	--	+-	+-

\*additional new data

The F.I.S.H. analyses gives an indication of binding sites for probes in the sandwich hybridization assay, but does not necessarily predict what combination of probes are best with respect to sensitivity and specificity. Previous experience has also shown that these combinations have to be determined empirically (Scholin C.A., Miller P.E., Buck K.R., Chavez, F.P., Harris, P., Haydock, P., Howard, J., Cangelosi G. 1998. Detection and quantification of *Pseudo-nitzschia australis* in cultured and natural populations using LSU rRNA-targeted probes. *Limnol. Oceanogr.* **42**: 1265-1272; Scholin unpublished work). The probes Het1.25F and Het 3F were chosen for the capture probes and the fluorescein labels were exchanged for a biotin label. An alternative probe, Het1.25bS, was also constructed to try and improve the specificity of Het1.25F which gave a weak signal when probed against other raphidophyte species in the F.I.S.H. format. Het1.25F is equivalent to Het1.25aS where the fluorescein label has been exchanged for biotin. The probes Het. sig2F, Raphid1F and Raphid2F were chosen for the signal probes and these probes were constructed with a fluorescein label.

### Example 7

#### Sandwich Hybridization Assays

Candidate probes were designed for the sandwich hybridization assay (S.H.A.) after screening the LSU rRNA gene with probes in a F.I.S.H. format as described in Example 6. Probes which displayed good specificity and sensitivity in the F.I.S.H. analyses were redesigned for the higher T<sub>m</sub> requirements of the sandwich hybridization assay, with biotin labels being attached for capture probes and fluorescein for the signal probes.

The sandwich hybridization assay utilized two probes. The first probe was a biotinylated species-specific capture probe, which was bound to a solid streptavidin-coated support. The capture probe then bound with target rRNA molecules and was transferred to a second solution containing a fluorescein-



labelled signal probe, thus forming the sandwich hybrid configuration. The hybrids were detected by a secondary-labelling reaction involving anti-fluorescein antibody conjugated to horseradish peroxidase that reacts with substrate to produce a blue colorimetric product. The intensity of the blue product is related to the number of rRNA molecules in the sandwich hybrid format, and from the number of rRNA molecules, it is possible to infer the number of algal cells present from a standard curve.

Samples for analysis with the S.H.A. were prepared as follows. Cultured cells and field samples were collected by gentle filtration onto 25 mm hydrophilic Durapore membranes (0.65  $\mu$ m pore size, Millipore). The filters were transferred to filter tubes (Porex, Fairburn, GA) containing 400  $\mu$ L of lysis buffer [(50 mM glycine, 10 mM EDTA, 5% N-lauryl sarcosine, 0.5% ProClin 150 (Rohm and Haas, Philadelphia), pH 11)], vortexed gently, then heated at 85°C for 5 min, with an additional vortex after 2.5 min of incubation. After the incubation, 600  $\mu$ L of hybridization buffer (100 mM Tris, 17 mM EDTA, 8.35% formamide, 5 M guanidine thiocyanate, pH 7.5) was added to tubes in heat block. The tubes were removed immediately, vortexed gently and allowed to cool for 5 min. The tubes were then capped with a filter tip (1  $\mu$ m pore size), or transferred to a syringe, coupled to a 0.45  $\mu$ m Durapore Millex filter unit (Millipore) to remove particulate material from the crude cell lysate. The cell lysates were either allowed to cool to room temperature (~5 min) then were processed immediately, or they were placed into a -80°C freezer for later analysis.

All hybridization steps were carried out using a portable robotics workstation which transfers a polystyrene strepavidin-coated strip with 12 prongs through each row of a standard 96-well microtiter plate (Saigene, Redmond, WA). Each row of the microtiter plate represents a different hybridization step. Assay development plates, which have all the reagents for the S.H.A., except for probes and samples, were provided by Saigene (Redmond, WA). Filtered cell lysate (200  $\mu$ L) was added to 3-4 wells for each sample, plus 2-3 wells for a negative control.

The negative control probes do not bind to any known rRNA molecule and serve as a procedural control to ensure that observed color development is specific to the capture/lysate/signal probe sandwich hybrid. The last well in the sample row was loaded with lysis/sample buffer. This well and the column serve as a chemistry-  
5 positive control, which utilises an oligonucleotide linker which forms a sandwich hybrid with the capture and signal probes.

After loading the samples into the wells, the microtiter plate was transferred to the temperature-controlled surface of the portable robotics workstation. The hybridization steps were carried out at 30°C and the entire run was completed in  
10 approximately 50 min. Each of the hybridization steps was completed with agitation of the prongs in an up and down motion in wells. Primary hybridization reactions, where the biotinylated capture probes are bounded to the strepavidin-coated prong, took place in 200 µL of Assay Wash buffer (Saigene, Redmond, WA), with capture probe at a final concentration of 200 ng µL<sup>-1</sup> for 10 min. The  
15 prong/capture probe hybrids were transferred to the sample wells and secondary hybridization was allowed to proceed for 10 min. After the secondary hybridization reaction, the prong/capture probe/analyte hybrids were transferred into wells containing 175 µL of Signal Buffer (Saigene, Redmond, WA), with signal probe at a final concentration of 300 ng µL<sup>-1</sup> for 8 min. At the completion of  
20 the tertiary hybridization step, the 'sandwich' hybrid (capture probe/analyte/signal probe) formation was complete. The sandwich hybrids were washed with 200 µL of Assay Wash solution (Saigene, Redmond, WA) for 2 min. Then the prong/sandwich hybrids were placed into 200 µL of anti-fluorescein-horseradish peroxidase (HRP) conjugate for 10 min. The prongs were then rinsed twice for 2  
25 min as described above. The final step was the incubation of prong/sandwich hybrids in 200 µL of HRP substrate (TMB) for 12 min to allow color development. At the completion of the run, the microtiter plate was quickly transferred into a plate reader and the absorbance was measured at 655 nm. The plate was removed and 50 µL of 10% (v/v) sulphuric acid was added for color enhancement. The plate

was transferred back into the microtiter plate reader and the absorbance was measured at 450 nm. All reagents for the S.H.A. are available from Saigene Corp. under the trade names listed above (Redmond, WA).

5 The sensitivity and cross-reactivity tests on cultured isolates for  
*Heterosigma akashiwo* revealed that the probe combination Het1.25aS/Raphid1F  
provided the best signal, and no cross-reactivity has been observed. It was  
expected that there may be some cross-reactivity with the Het1.25aS probe based  
on the observations of the probe Het1.25aF in the F.I.S.H. analyses. The higher  
stringency and the lengthening of the probe for the S.H.A. appears to have  
10 eliminated these potential problems. It is well known that a small shift in a probe's  
position can have dramatic effects on the signal output. The probe Het1.25bS  
which was shifted a few bases towards the 5' end of the LSU rRNA gene, and  
which still included numerous bases of the probe Het1.25aS had a significantly  
reduced signal output in the S.H.A. in comparison to the Het1.25aS/Raphid1F  
15 probe combination. The results for the sandwich hybridization assays for the probe  
combination matrix are described in Table 5.

**Table 5. Sandwich hybridization probe matrix showing the combination and position of probes, and the signal output at 450 and 655nm.**

Probe Combination	Target	Probe Positions	OD@450nm	OD@655nm
	Species <sup>b</sup>	5' end <sup>a</sup>	±SD	±SD
Het1.25aS/Raphid1F	<i>H. akashiwo</i>	149 / 79	1.595 ±0.209	0.768 ±0.091
Het1.25aS/Raphid2F	<i>H. akashiwo</i>	149 / 344	0.0633 ±0.006	0.058 ±0.003
Het1.25aS/Het.sig2-3'F	<i>H. akashiwo</i>	149 / 366	0.0773 ±0.006	0.063 ±0.002
Het1.25bS/Raphid1F	<i>H. akashiwo</i>	157 / 79	0.583 ±0.081	0.307 ±0.040
✗ Het1.25bS/Raphid2F	<i>H. akashiwo</i>	157 / 344	0.094 ±0.048	0.073 ±0.024
Het1.25bS/Het.sig2-3'F	<i>H. akashiwo</i>	157 / 366	0.054 ±0.003	0.054 ±0.002
Het3S/Raphid1F	<i>H. akashiwo</i>	567 / 79	0.093 ±0.009	0.068 ±0.003
Het3S/Raphid2F	<i>H. akashiwo</i>	567 / 344	0.222 ±0.011	0.128 ±0.002
Het3S/Het.sig2-3'F	<i>H. akashiwo</i>	567 / 366	0.080 ±0.016	0.080 ±0.016
Fib1.25aS/Raphid1F	<i>F. japonica</i>	122 / 49	1.32 ±0.040	0.634 ±0.020
Fib1.25aS/Raphid2F	<i>F. japonica</i>	122 / 322	0.094 ±0.009	0.067 ±0.004
Fib1.25aS/Fib.sig2-3'F	<i>F. japonica</i>	122 / 344	0.087 ±0.006	0.067 ±0.004

<sup>a</sup>Ribosomal RNA complement. Aligned position with respect to target sequence.

<sup>b</sup>The clones used for this analysis are CAWR05 FOR *Heterosigma akashiwo* and CAWRO3 for *Fibrocapsa japonica*. Refer to Table 1 for details about these clones.

5

The best combinations of probes for signal production in *Heterosigma akashiwo* and *Fibrocapsa japonica* were the Het1.25a capture probe/Raphid1F signal probe and Fib1.25a capture probe/Raphid1F signal probe respectively. These combinations of probes were used in all the following analyses.

5           Also, the S.H.A. probes were modified from the F.I.S.H. probes due to the higher Tm requirement and this may affect the signal output and specificity. A probe matrix was used to test all possible combinations of capture and signal probes for signal production. The best combination of probes for *Heterosigma akashiwo* and *Fibrocapsa japonica* were then tested for specificity in the S.H.A.  
10           format. Figure 1 shows specificity and isolate comparison for the sandwich hybridization assay combinations targeted at *H. akashiwo* and *F. japonica*. A). The probe combination Het1.25aS/Raphid1F which is targeted at *H. akashiwo* was screened against a number of *H. akashiwo* isolates and other raphidophyte species. *H. akashiwo* 1 represents the isolate *H. akashiwo* CAWR04; *H. akashiwo* 2 =  
15           CAWR05; *H. akashiwo* 3 = CAWR09 and *H. akashiwo* 4 = CAWR14. The *F. japonica* isolate was CAWR02. B). The probe combination Fib1.25aS/Raphid1F which is targetted at *F. japonica* was screened against two *F. japonica* isolates and other raphidophyte species. *F. japonica* 1 represents the isolate CAWR02 and *F. japonica* 2 represents the isolate CAWR03. The  
20           *H. akashiwo* isolate was CAWR05.

#### Example 8

##### Standard Curve Construction

25           The standard curves were constructed by spiking cells which were in log phase growth into 1 L of either filtered sea water (FSW) or field sample, then collecting the cells by gentle filtration and lysing as described above. The concentrated stock solutions were then serially diluted with either lysis/hybridization buffer for the FSW sample or field background for the field

sample. Four replicates of each dilution were loaded onto the hybridization plates to check for variation of signal from well to well.

Figure 2 shows serial dilution standard curves for *F. japonica* and *H. akashiwo* based on the sandwich hybridization assay; A) standard curve at 450 nm for *F. japonica* cells spiked into field and filtered sea water (FSW) background; B) standard curve at 450nm for *H. akashiwo* cells spiked into field and FSW background.

The differences between the spiked-FSW and field sample standard curves for both *Heterosigma akashiwo* and *Fibrocapsa japonica*, were non-significant, indicating that a complex field background has minimal impact with regard to efficacy and sensitivity of the S.H.A.'s for *Heterosigma akashiwo* and *Fibrocapsa japonica*.

#### Example 9

##### Signal Comparison Over a Complete Growth Cycle

The signal generated by cells through a growth cycle was examined by innoculating two batch cultures per species (1 L of f/2 media in 5 L flasks) with 1000 *Heterosigma akashiwo* cells mL<sup>-1</sup> and 500 *Fibrocapsa japonica* cells mL<sup>-1</sup> respectively. These cultures were grown under a 12:12 LD photoperiod at 20±1°C. Samples were collected at regular intervals covering the entire growth cycle. The *Heterosigma akashiwo* samples were preserved with Lugols Iodine, and at least three 10 µL aliquots were counted. *Fibrocapsa japonica* was preserved with the ethanol/saline fixative and counted as for *Heterosigma akashiwo*. Lugols Iodine preserved *Fibrocapsa japonica* resulted in clumped cells which were not quantifiable. The S.H.A. lysates were prepared by collecting 20,000 *Heterosigma akashiwo* cells 1 ml lysate<sup>-1</sup> and 10,000 *Fibrocapsa japonica* cells 1 ml lysate<sup>-1</sup>, which corresponds to 4000 and 2000 cells well<sup>-1</sup> (well = 0.2 ml) respectively.

These samples were stored at -80°C and once all samples were collected they were screened with the S.H.A. The samples were processed in one lot to minimize any differences in various batches of S.H.A. plates and reagents.

5 There was no significant difference between the FSW and field treatment standard curves for either *Heterosigma akashiwo* or *Fibrocapsa japonica*. The results for *Heterosigma akashiwo* and *Fibrocapsa japonica* are displayed in Figure 1.

10 Preferred long-term storage of lysate is considered to be storage of samples at -80°C for extended periods of time. A large batch of lysate was prepared, split into aliquots and frozen at -80°C to check the effect of storage on signal. At various times over a year, a sample was removed and analyzed to check for degradation and subsequent signal deduction and the effect of multiple freeze/thaw cycles was tested although data are not provided here.

15 Samples which were required for later analyses were stored at -80°C. These samples were immediately stored either as a lysate or as cells filtered down onto a membrane. The lysates showed no appreciable drop in signal over an entire year of storage, whereas the processing of cells which had been filtered onto membranes resulted in a substantial decrease in signal. Repeated freeze/thaw cycles of the lysate also resulted in the rapid decrease of signal.

20 The growth cycle experiment showed a 2-fold signal variation for *Heterosigma akashiwo* between the two batches in the first 9 days of the culture, but after day 9 the signal from the batches was comparable. In late stationary phase where the cultures were starting to die, there was a 2-fold decrease in signal. For *Fibrocapsa japonica* there was a variance of about 2-fold over the entire growth  
25 cycle with no appreciable dropoff in signal at late stationary phase growth. The results for *Heterosigma akashiwo* and for *Fibrocapsa japonica* are displayed in Figure 2.

Figure 3 shows ribosomal RNA levels as measured using the sandwich hybridization assay over an entire growth cycle for *H. akashiwo* and *F. japonica*. In the graphs B and D the open box data points represent signal at 450 nm, whereas the closed circles represent signal at 655 nm; A) Changes in cell density over a growth cycle for a batch culture of *H. akashiwo*; B) Signal from the sandwich hybridization assay for a constant number of cells over a growth cycle of *H. akashiwo*; C) Changes in cell density over a growth cycle for a batch culture of *F. japonica*; D) Signal from the sandwich hybridization assay for a constant number of cells over a growth cycle of *F. japonica*.

The determinations of signal from cells grown over an entire growth cycle using batch cultures showed that the signal only varied by a factor of 2, except for dying cells in late stationary phase growth for *Heterosigma akashiwo*. Using F.I.S.H. analyses there is an order of magnitude decrease in signal for cells entering late stationary phase growth. The difference in observations for F.I.S.H. and S.H.A. formats suggests that the rRNA in preserved cells may not be as accessible due to the increase in polysaccharides in cells and/or the protein interaction with rRNA. Further research into these observations is required to gain insight into the fate of rRNA pool during a growth cycle. Further, these observations are based on cultured lysates and may not reflect the situation in the field. The variation between signal in the *Heterosigma akashiwo* batch cultures is evident in the first 9 days of growth. These differences may be due to inaccurate counts and it was also observed that one of the cultures had significant aggregation/sticking of cells to side of flask during these 9 days.

No non-specific signal production for either combination of probes has been observed for cultured species or field samples. The probe combinations for *Heterosigma akashiwo* and *Fibrocapsa japonica* gave a signal from the S.H.A. which was comparable for all geographic isolates of the target species.



## Example 11

### Field Trials

The preliminary field trials using the S.H.A. have detected *Fibrocapsa japonica* in low numbers (2-3000 cells L<sup>-1</sup>, which corresponds to 160 and 240 cells well<sup>-1</sup> respectively) at three locations in the Hauraki Gulf, East Coast of the North Island, New Zealand. These results were confirmed independently by the commercial monitoring team of the Cawthron Institute, Nelson, New Zealand, using traditional methodology. One unconfirmed 'positive' result has been registered for *Heterosigma akashiwo* in a sample collected from the Santa Cruz pier, Santa Cruz, California.

Field trials were conducted on sandwich hybridization assays for *Heterosigma akashiwo* and *Fibrocapsa japonica* using opportunistic samples sent to the Cawthron Institute commercial monitoring team, Nelson, New Zealand. The sandwich hybridization assay for *Heterosigma akashiwo* and *Fibrocapsa japonica* was compared against standard light microscopy counts. Their counts were based on a single 10 mL sub-sample of the acid Lugol's Iodine preserved samples. Overall, comparisons of cell numbers as counted by light microscopy and the S.H.A. are in good agreement. No false positives were observed and all samples that contained either *Heterosigma akashiwo* or *Fibrocapsa japonica* above the lower limit of detection were detected using the S.H.A. This observation indicates that New Zealand populations of these two species have minimal genetic diversity at the Large Subunit rRNA (LSU rRNA) gene level. Also, work on cultured isolates from Australia, Canada, Europe, Japan, Korea, and the USA indicate that global populations of *Heterosigma akashiwo* and *Fibrocapsa japonica* are relatively homologous based on the LSU rRNA gene. This result is extremely beneficial for the deployment of this assay for these species globally. In a recent paper submitted for the proceedings of the 9th International Conference on Harmful Algal Blooms Conference, Tasmania, Australia, Rhodes et al. Concluded

PATENT

Our Ref: 506812000120

“...it is clear that this assay is now ready to be integrated into a suite of monitoring tools for both fish farm managers ...”.

- 43 -